

Specification

Filtration · Rapid Tests · Water Analysis · Chromatography · Bioanalysis
Filtration · Schnellteste · Wasseranalytik · Chromatographie · Bioanalytik


Product:	CHROMAFIL [®] Xtra Polyamide (PA) / Nylon
REF:	729213
	Filters are HPLC certified for low level UV-absorbing substances for gaining accurate chromatographic results
Housing material:	Polypropylene (PP)
Labeling:	membrane type imprinted on filter top
Membrane filter material:	Polyamide (Nylon)
Pore diameter:	0.45 µm
Filter diameter:	25 mm
Filtration area:	~ 4.9 cm ²
Pressure stability:	> 12 bar
Connections inlet:	Female luer lock
Connections outlet:	Male luer
Dead volume:	~ 80 µL

HPLC-test


Conditions:	Two mL of the solvent (specified under the chromatograms) were applied to the filter; 100 µL of the filtrate were injected to the HPLC.
Eluent A:	water
Eluent B:	acetonitrile
Gradient:	10% → 95% B in 25 minutes
Flow rate:	0.5 mL/min
Sensitivity:	-5 to 10 mAU at 254 nm
Column:	125/4 NUCLEODUR [®] C18 Gravity 5 µm (REF: 760100.40)

Specification


Filtration · Rapid Tests · Water Analysis · Chromatography · Bioanalysis
Filtration · Schnellteste · Wasseranalytik · Chromatographie · Bioanalytik

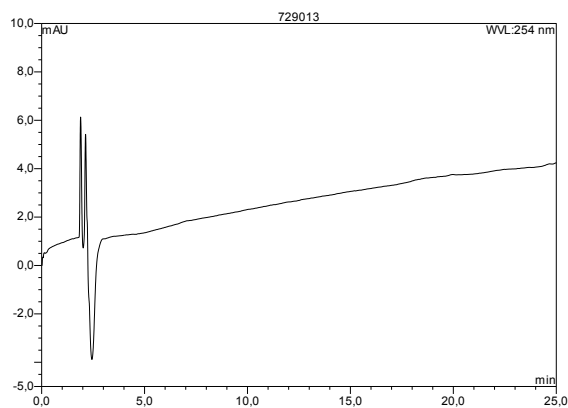

Methanol filtrate

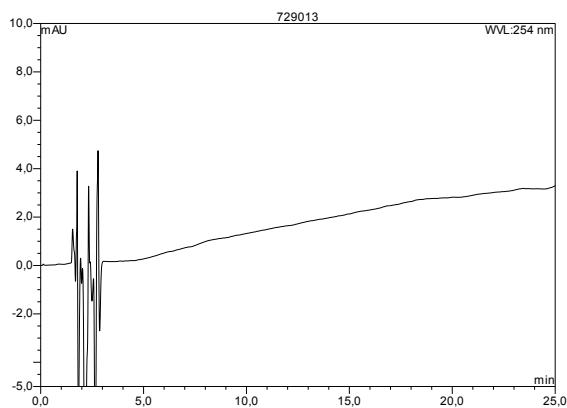

Acetonitrile filtrate

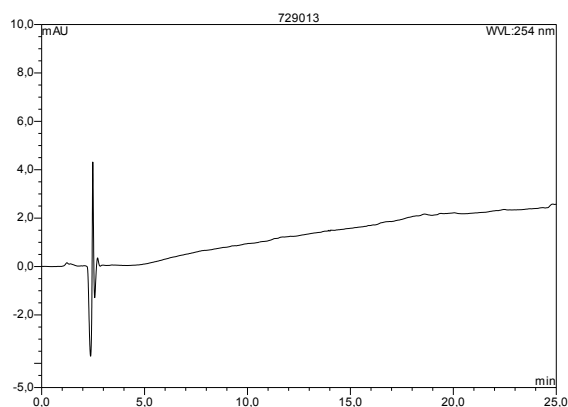

Water filtrate

There are no impurities found > 0.0001 % as potassium biphthalate.

